




ORIGIN[®] DUOTRAN[™] PRO

www.origincadcam.com | 800-255-2839 | 801-281-4012 | 2371 So. Presidents Dr., Ste. E, West Valley City, UT 84120

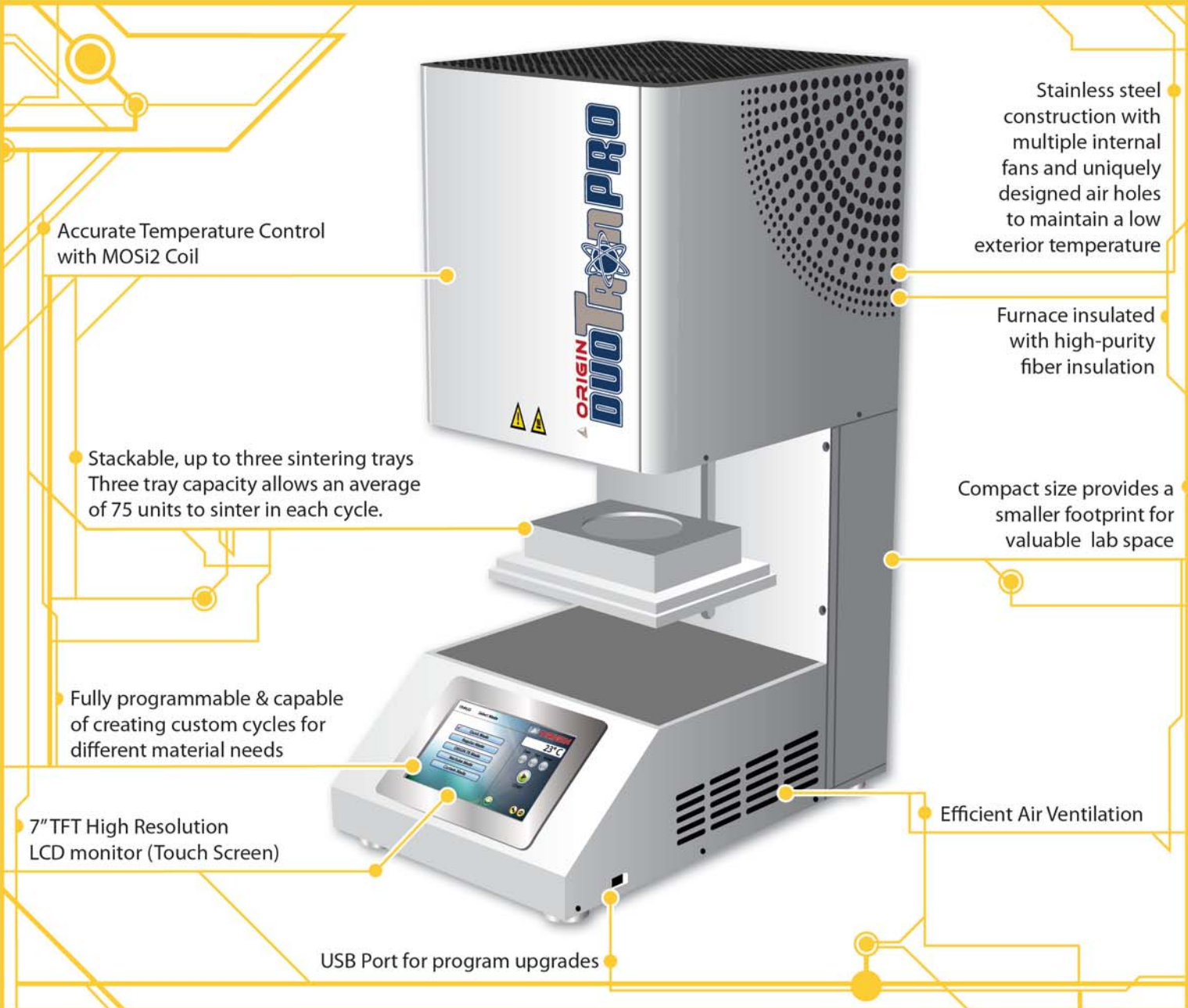

Accurate Temperature Control with MOSi2 Coil

Stainless steel construction with multiple internal fans and uniquely designed air holes to maintain a low exterior temperature

Furnace insulated with high-purity fiber insulation

Compact size provides a smaller footprint for valuable lab space

Stackable, up to three sintering trays
Three tray capacity allows an average of 75 units to sinter in each cycle.

Fully programmable & capable of creating custom cycles for different material needs

Efficient Air Ventilation

7" TFT High Resolution LCD monitor (Touch Screen)

USB Port for program upgrades

Fast Sintering!

\$16,900⁰⁰

- Increases zirconia production speeds & decreases production costs
- Cost Savings and Environmental Friendliness
- Compact, requiring less equipment space
- Consistent zirconia color
- Eliminate contamination build up

Mode & Indication

Quick Mode: Single copings or small bridges (3 to 4 unit bridges)

Regular Mode: Large bridge cases or 4 unit bridge cases with large pontics (4 - 7 hours)

Full Contour Mode: Full contour zirconia sintering at an increased temperature.

Custom Mode: Build your own sintering cycle to fulfill all your zirconia needs (LAVA, Cercon, Prismatic, etc)